



Kiss

Veneering Porcelain


# Craft meets creativity

Brochure for the dental laboratory



For all framework materials

# 3 Kiss ceramic lines

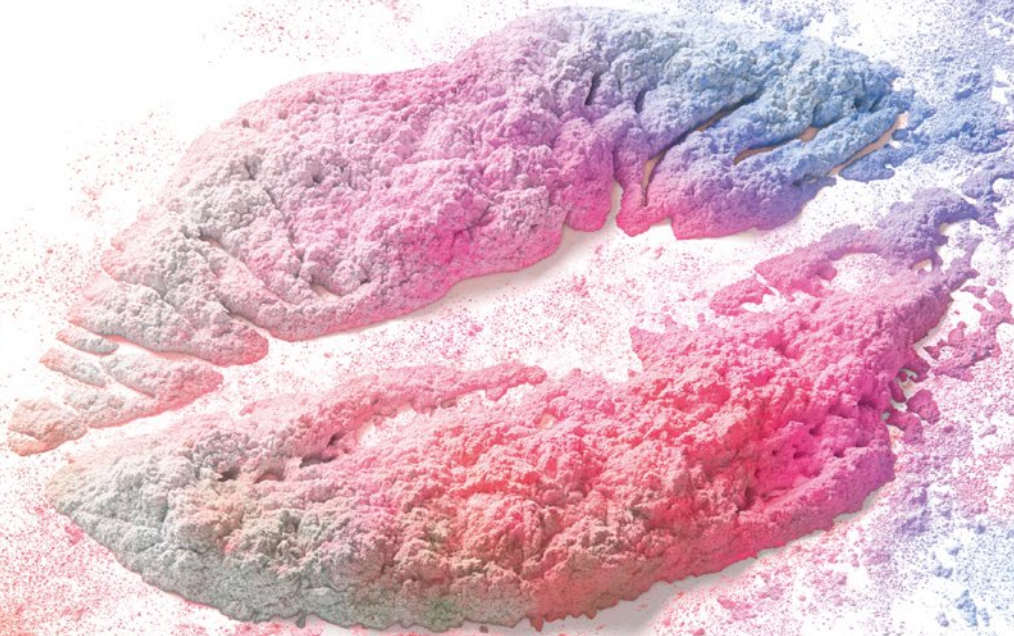

---

## DUCERAGOLD KISS

High-expansion alloys

---

## DUCERAM KISS

Classic alloys, bio-alloys,  
non-precious alloys

---

## CERCON CERAM KISS

Zirconia (Cercon base, Cercon ht)

- One build-up scheme for all framework materials
- Simple and fast processing thanks to the perfect adaptation of the ceramic veneering material to the system's framework materials
- Each ceramic veneering material within the Kiss system is associated with closely matched framework materials offering you a maximum of processing safety
- Excellent shade harmony across the Kiss ceramics
- Multi-functional and light-dynamic materials powders
- Efficient ceramic range – only 73 different materials powders

Keep it simple and safe

# Kiss veneering ceramics



## Reduced number of powders

### + **Reduced material concept without sacrificing quality**

The accurate categorization of the initial and intensive powders, the tinted opaque and dentins as well as the universal mixing scheme make it possible to achieve highly aesthetic veneers without any quality loss: Use only 73 materials to do 100% of all jobs.

### + **Economic advantage**

The reduced number of materials benefits your bottom line.

## Perfect aesthetics

### + **More shade fidelity**

Optimized shade pigments enhance the shade fidelity of Kiss, similar to the V-shades. This also largely eliminates the problem of differences in shade appearance under different lighting conditions.

### + **Opalescent effect remains intact**

In the classic high-fusing veneering ceramics, the opalescence can lose its effect after a number of firing cycles. With Kiss, this effect – and with it the natural light dynamic appearance – remains intact thanks to the patented manufacturing process used for the Kiss opal powders.

## Simple and safe

### + **Simple and safe shade production and reproduction**

Secure processing using a simple layering system for basic and customized build-ups.

### + **Perfectly matched shades across ceramic lines**

Easy and problem-free processing of different framework materials possible in the same patient case.

### + **Uniform build-up system for all framework materials**

No extra training cost or learning curve when switching from one framework material to another.

The Kiss veneering concept

# Materials and shades

**Economical 3-mass build-up**  
Basic build-up for V-shades

**Intelligent mixing for highly varied shades**  
Power Chroma share recognition

TC Opaker Orange

## 100% of the cases with only 40 materials

All V-shades can be produced with the 40 materials shown here, without any mixing at all. The range of six incisals allows a better representation of the V-shades in the incisal area. All opaques are precisely adjusted to the respective basic shade. Their fluorescence levels are adapted to the different shade intensities and, consequently, to the natural tooth model.


A superb base for your daily work – effortless success for almost everything you attempt.

## Fluorescent Power Chromas


The six fluorescent Power Chromas can reproduce the majority of all cervical and mamelon effects as well as increases in chroma. Simple 1:1 mixtures between these Power Chromas results in an additional 15 intermediate shades.

## Natural opalescence and fluorescence effects

### Incisals for individual value control

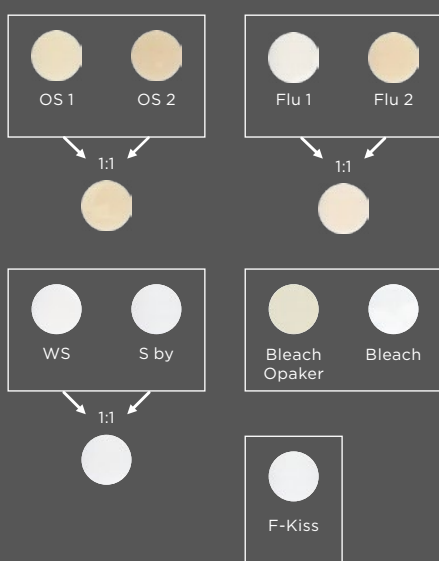

### Opal effects + multifunctional Stand by materials

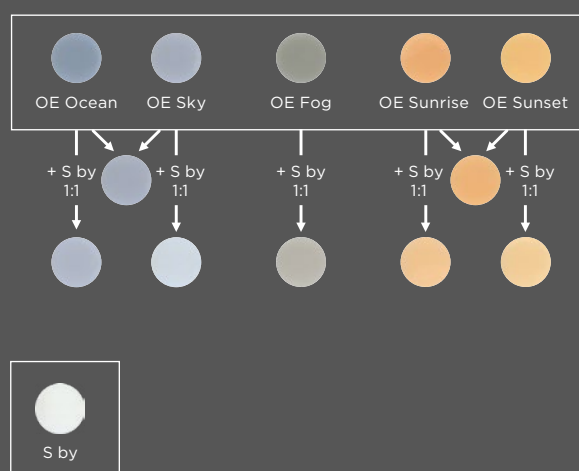

Depicted colours may vary from the original due to the printing method.

## Custom incisals, opal effects and multifunctional Stand by powders

The individual incisals within the Kiss concept allow the emulation of all natural opalescent and fluorescent effects. Opal incisals 1 and 2 are intended for lighter (OS 1) and darker (OS 2) shades. For intermediate shades, the two powders are simply mixed at a 1:1 ratio. The same principle holds for Flu Inside 1 and 2. These highly fluorescent materials are designed for the inner layers, covering the opaque even where space is limited, while at the same time increasing the lightness in the incisal area. White Surface (WS) is a whitish opalescent effect used to increase surface value (brightness). The effect can be minimized

by adding Stand by. To reproduce the effects of teeth that have been bleached, Bleach Opaques and Bleach Dentins are available. A particularly successful and versatile correction material is the transparent Final Kiss. The Stand by multifunctional material is a strongly transparent opal effect that is one of the keystones of the Kiss concept. It can be used by itself or mixed with any other material. The three opal effects, Ocean, Sky and Fog, control value and opalescence in bluish and greyish incisal regions. Sunset and Sunrise are suitable for incisal characterization and for subdued transparent chroma gradations in the body region.




The Kiss veneering concept

# Materials and shades

## 4 Gum materials

for individual gingival tints

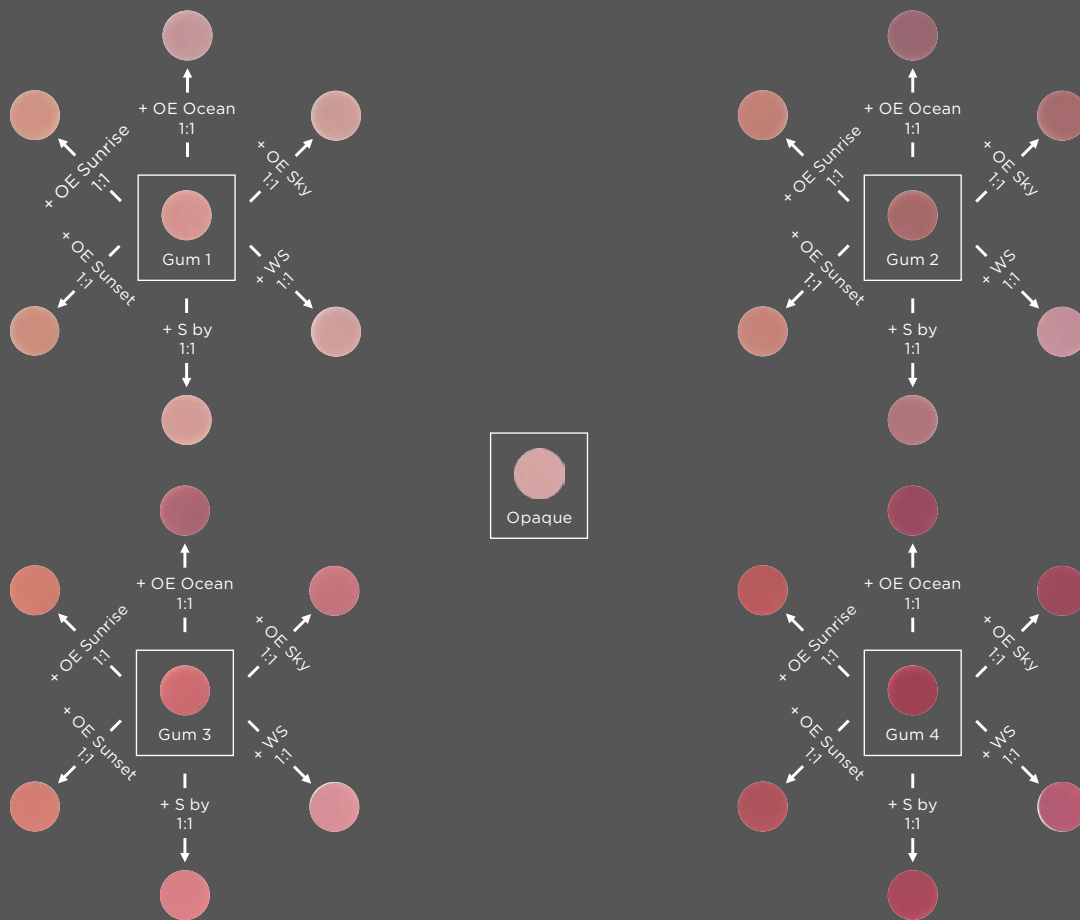


## 4 gingival materials

By mixing Gum 1 to Gum 4 with selected effects (see chart), natural gingival effects can be emulated very closely and simply.


Basic or custom build-up technique

# Aesthetics without limitations

Basic build-up  
in the anterior region

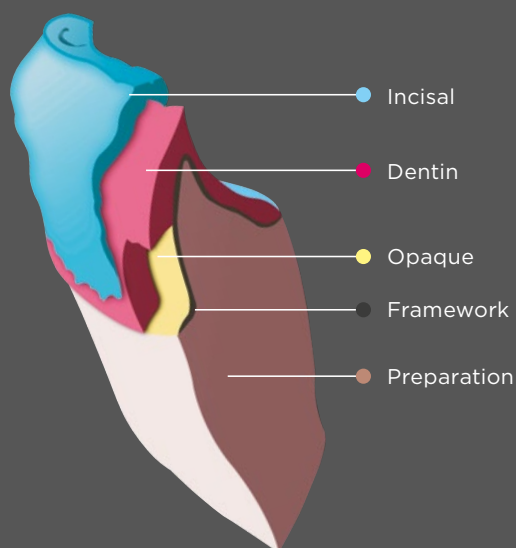

Custom build-up technique  
in the anterior region

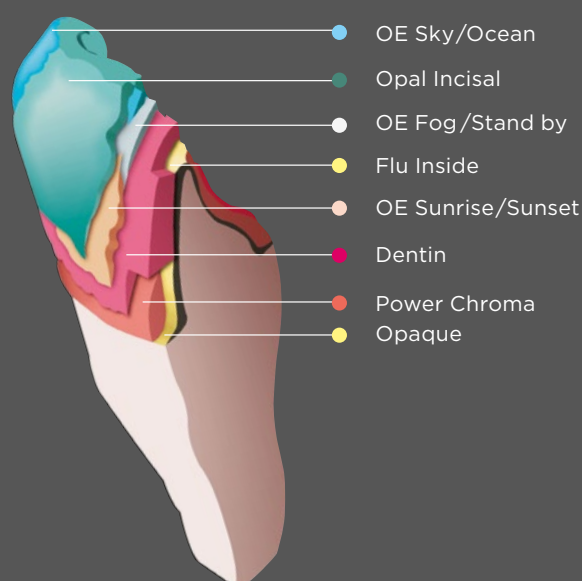

- Use of Power Chroma in the cervical area for customization
- Dentin build-up



- Flu Inside in the incisal third for value control
- An incisal-edge overlay with OE Ocean or Sky for a youthful, bluish incisal ridge



- In the body area OE, Sunset and Sunrise
- In the incisal region, opal incisors 1 or 2 and as well as Transpa

The ceramic line

# Cercon ceram Kiss

Aesthetic/cosmetic dentistry is intricately linked with all-ceramic dental restorations. To meet expectations in this high-end segment, Cercon ceram Kiss has been specifically developed to adapt to zirconia with a CTE of  $10.5 \mu\text{m}/\text{m} \cdot \text{K}$  (25-500 °C).

Cercon ceram Kiss is equally well suited for use on Cercon base and Cercon ht frameworks. Cercon ceram Kiss lets you create highly aesthetic ceramic veneers, simply and quickly.





---

## CERCON CERAM KISS FEATURES

- + Can be processed with or without a liner
- + Homogeneous firing results
- + Opal incisals for maximum aesthetics
- + Safety proven through chewing simulator tests
- + Easy handling thanks to our familiar Kiss materials and build-up philosophy



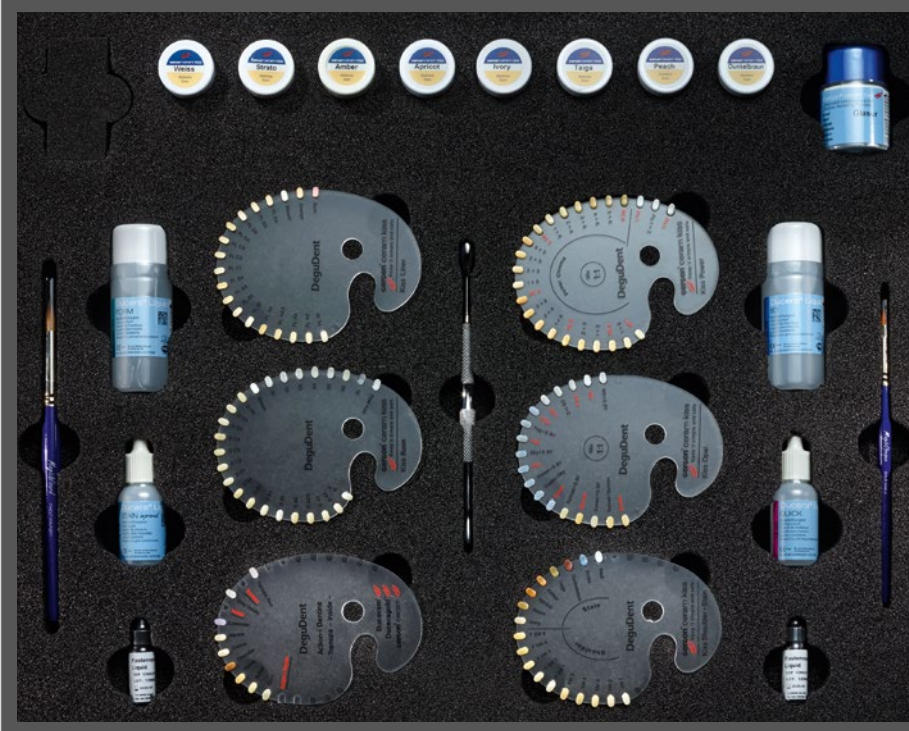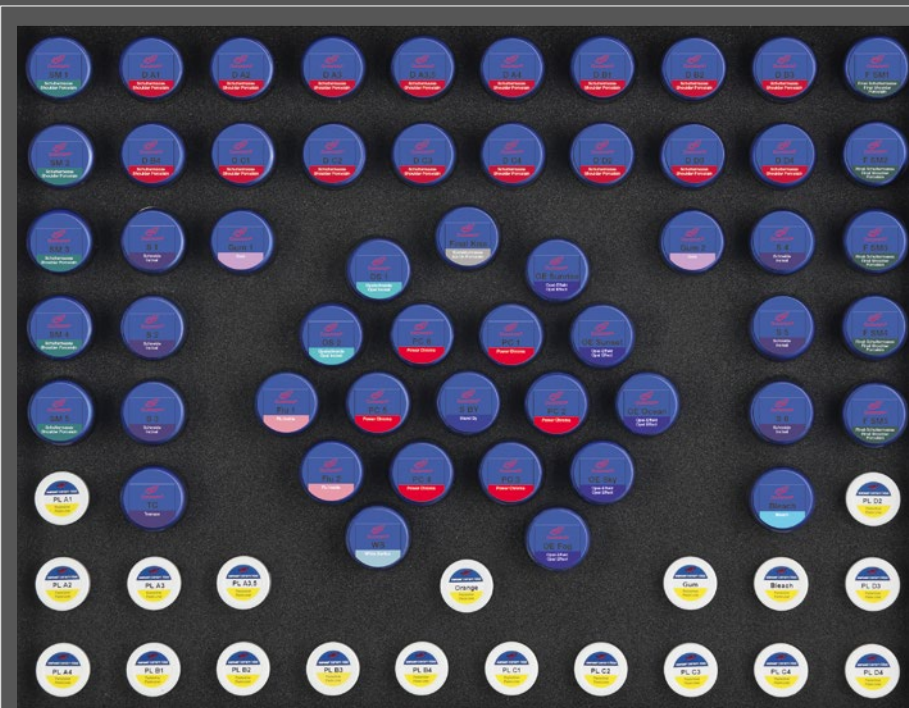

The ceramic line

# Duceragold Kiss

Duceragold Kiss is the ceramic veneering material for high-expansion alloys with a CTE of  $16.2$  to  $17.6 \mu\text{m}/\text{m} \cdot \text{K}$  ( $25$ - $600^\circ\text{C}$ ). It is a component of the GoldenGate System, widely known and appreciated for more than 20 years.

Duceragold Kiss is also indicated for veneering Cergo Kiss pressable ceramics. A dual-benefit veneering concept!




---

## DUCERAGOLD KISS FEATURES

- + Antagonist-friendly
- + Homogenous surfaces
- + Shade variety through simple 1 : 1 mixtures
- + Light dynamic powders
- + For all-ceramic and metal-ceramic restorations



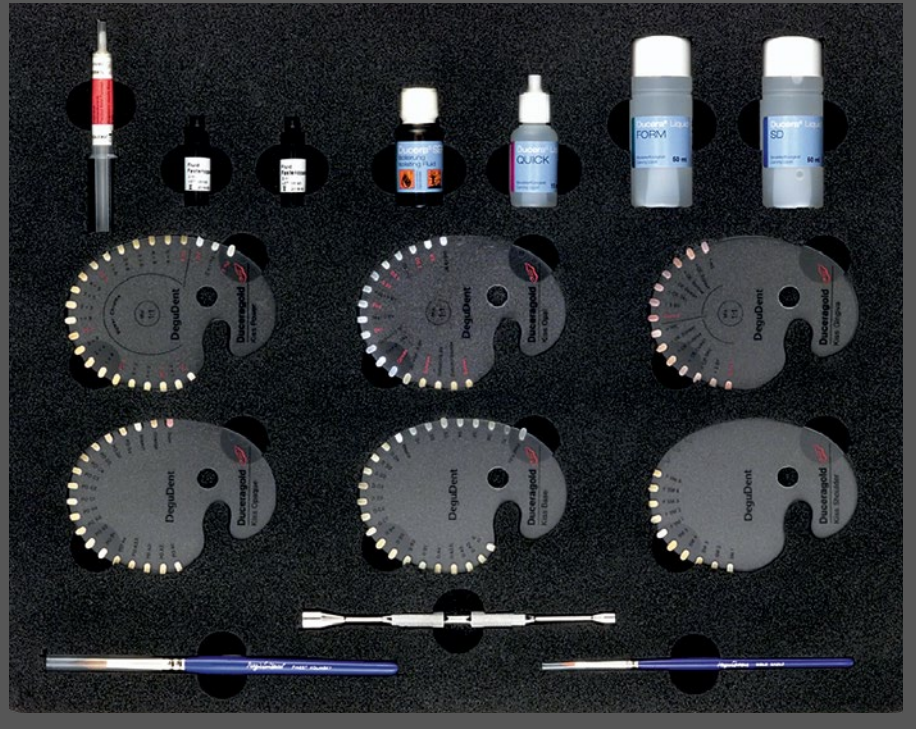



The ceramic line

# Duceram Kiss

Duceram Kiss is the ceramic veneering material for classic precious and non-precious alloy frameworks with a CTE of  $13.8$  to  $15.4 \mu\text{m}/\text{m} \cdot \text{K}$  (25-600 °C).

Duceram Kiss thus covers a broad range of indications, starting from the copper- and palladium-free alloys, BiOcclus Kiss, and classical framework alloys including non-precious alloys such as StarLoy soft and Duceralloy C+S.



The most distinctive characteristics of Duceram Kiss are its excellent shade fidelity, robust and simple handling and firing stable opal incisals, ensuring maximum aesthetics.



---

## DUCERAM KISS FEATURES

- + Firing stable Opal incisals for maximum aesthetics
- + Simple and accurate reproduction of V shades
- + Three powder standard build-up technique for maximum economy
- + Maximum shade variety

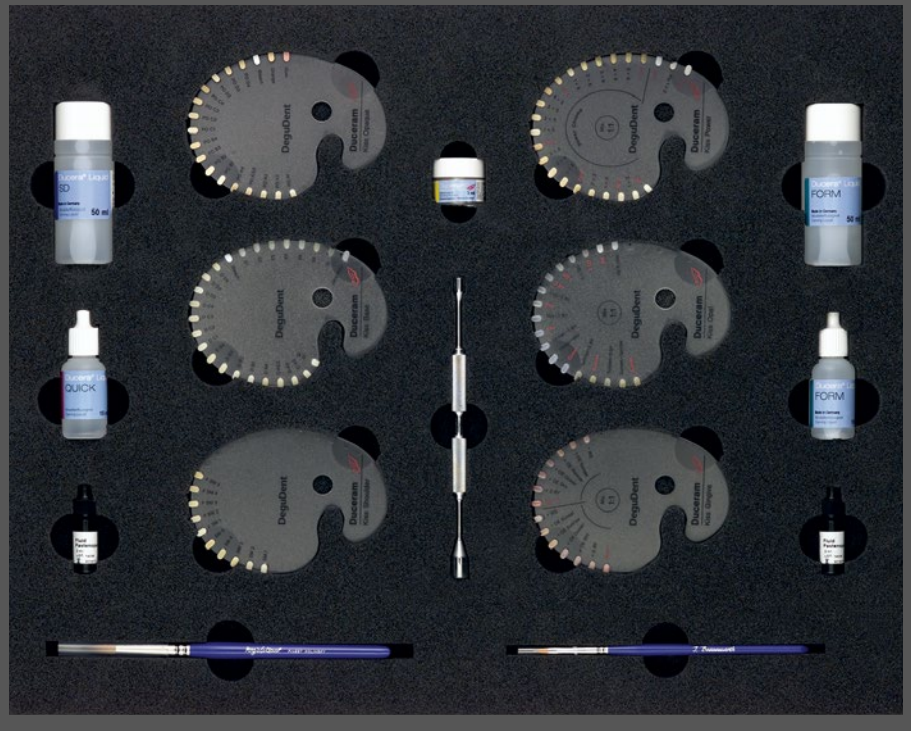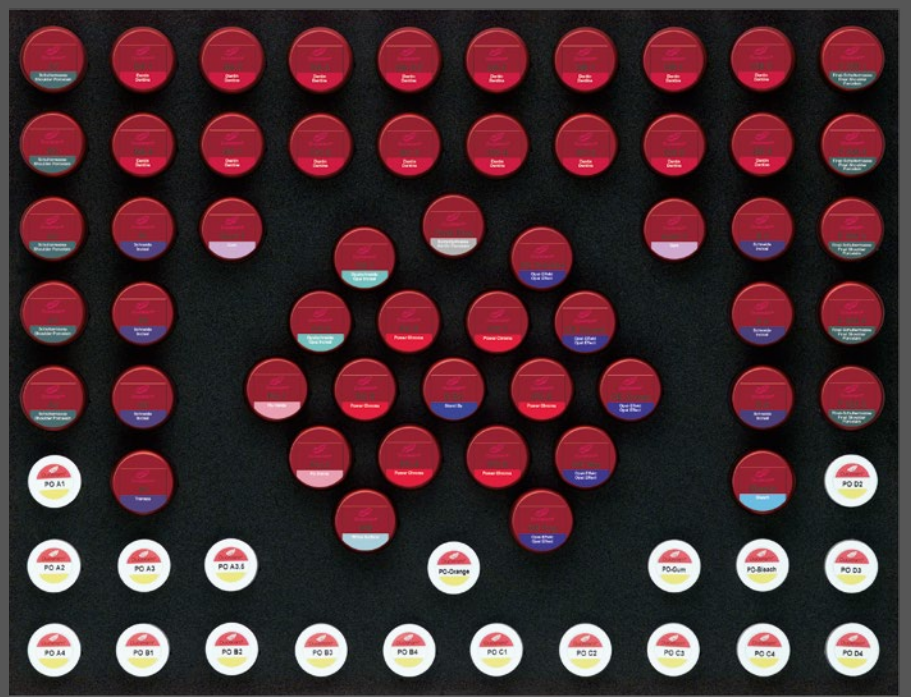

A complete system developed to perfection

# New: Kiss Artist Kit

## Action-i Dentine



### Corn

Translucent; similar to opal incisals to support D-shades, also suitable for alternating layers for mamelons.



### Butter

Translucent; similar to opal incisals to support B-shades, also suitable for alternating layers for mamelons.



### Honey

Translucent; similar to opal incisals to support A-shades, also suitable for alternating layers for mamelons.



### Lavender

Translucent; similar to opal incisals to support C-shades, also suitable for alternating layers for mamelons.



### Creme

Modifier; to be applied in the upper third of the dentine layer. Can be used pure or mixed.



### Marble

Modifier; to be applied in the upper third of the dentine layer. Can be used pure or mixed.



### Chocolat

Modifier for additional darkening of darker tooth shades; to be applied to the entire body, can be used pure or mixed.

## Fluorescent Dentine



### Gray Inside (GI)

A unique fluorescent dentine that can be used to compensate for high-value framework contours in the incisal region.

## Transpa effects



### Transpa Red (TR)

For additional customization of reddish tooth shades.

## Gums



### Gum 3

Complement to Gum 1 and Gum 2 - chromatic.



### Gum 4

Complement to Gum 1 and Gum 2 - highly chromatic.



# Kiss Action-i Dentine





DeguDent GmbH  
Rodenbacher Chaussee 4  
63457 Hanau-Wolfgang  
Germany  
[www.dentsplysirona.com](http://www.dentsplysirona.com)

22269/REV 2017-05

