


THE DENTAL
SOLUTIONS
COMPANY™



Teneo

Innovation in perfect harmony

dentsplysirona.com/teneo

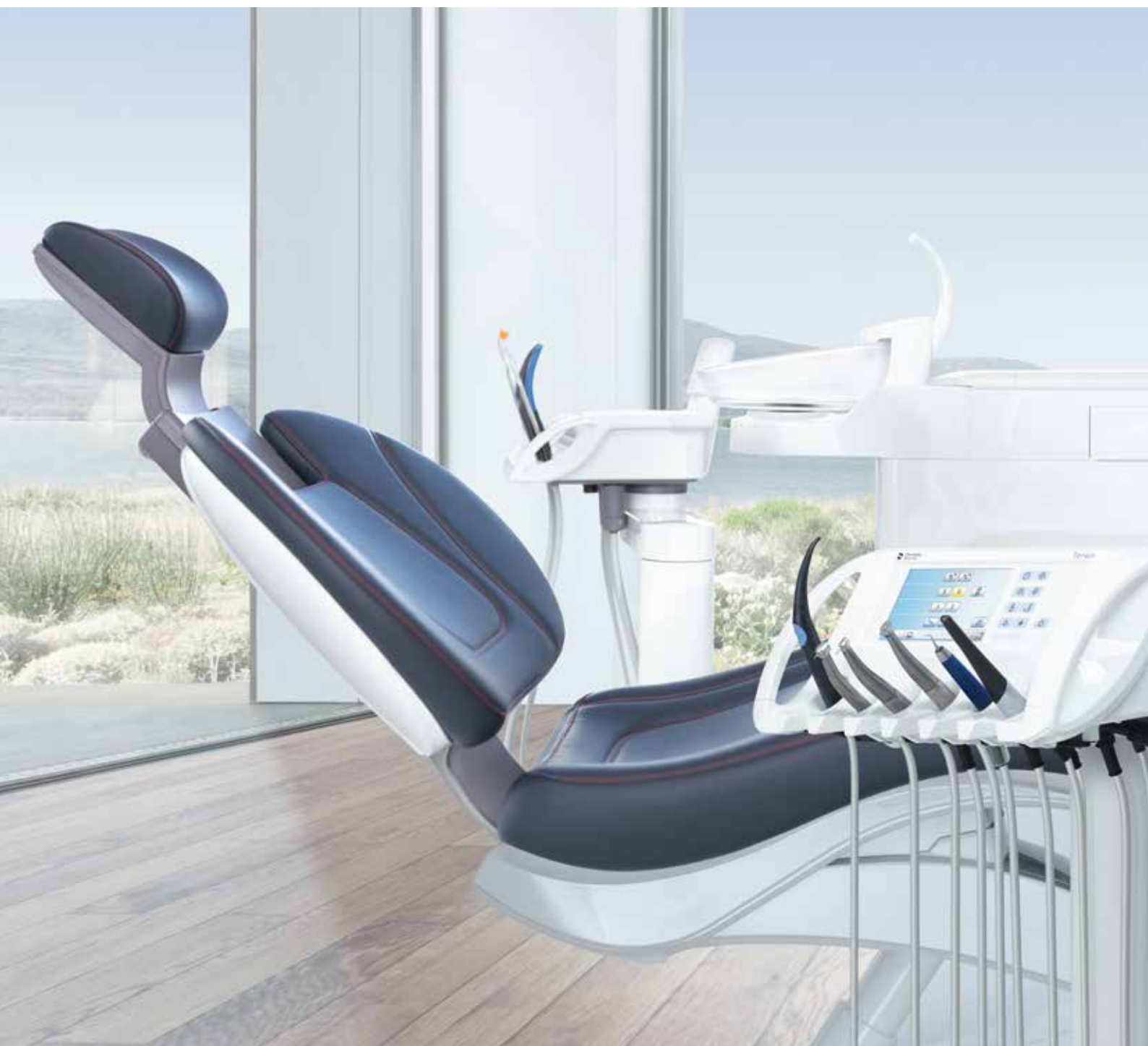



Innovation in perfect harmony

Teneo takes an uncompromising approach: Only the best is good enough. You can achieve a perfectly harmonious workflow if you and your treatment center interact optimally. For this reason, Teneo responds automatically and individually to all your needs. Carefully selected features assist you in every treatment situation, down to the most minute detail. Comfort features ensure that your patients feel more relaxed throughout the treatment. So you can focus entirely on your work. Why settle for less?



Treatment redefined

Your daily treatment efforts are the key to the success of your practice. Teneo has been developed to consistently match the dental workflow. Only if you can work efficiently and ergonomically your patients will benefit from the full range of your professional abilities. Teneo provides automatic and motorized options to adapt completely to your needs. Teneo responds directly to you and your patients – and the next treatment step is only one small movement away.

Why settle for less?

With motorized features such as the motorized sliding track, motorized headrest, and motorized cuspidor, Teneo adjusts perfectly to your workflow.

You know every move – so does Teneo

To let you focus fully on your patients, Teneo and its motorized sliding track places everything you need at your fingertips – hands-free and effortlessly. You and your assistant can maintain an ergonomic working position at any time.



Split concept for a better handling range



The dentist and assistant elements as well as the tray can be positioned independently to let you switch flexibly between team and solo treatment modes.

Hands-free operation



The wireless foot switch is your multifunctional control element for all chair functions. It provides hands-free control of the EasyTouch display and the motorized sliding track.

You want to treat ergonomically – Teneo adapts

A healthy operator posture during treatment requires the patient to be correctly positioned and maximum access ensured. Teneo provides comprehensive treatment comfort for dentists and patients alike, with its motorized headrest, cuspidor, and sliding track.



Vertical mouth point positioning



Teneo knows your ideal working posture. So that you can always work in your optimum position, the vertical mouth point is always positioned at the same working height — even with patients of different stature (see illustrations).

Motorized headrest



For optimum access to the treatment site, the extension and inclination of the motorized headrest can be adjusted using the foot switch, the user interface, or even manually.

Motorized cuspidor



The cuspidor is motorized and swivels in and out. Its position is programmable, contributing to fast and safe procedures.

The difference is in the details

Most of what any treatment is about is to keep the flow intact. For that you need a treatment center whose functions are meticulously adapted to your workflow — and not the other way around. So Teneo comes to your treatment room with an extensive range of functions already implemented. The user interface is fully customizable to let you apply the full range of treatment modalities.

Why settle for less?

Teneo is ideally preconfigured with a generous user interface, 6 + 1 customizable instrument positions and numerous aids for your patient communication. With one keystroke it is perfectly suitable for endodontic or implant procedures.




You want total freedom – Teneo is flexible

Integrated treatment functions for implantology and endodontics make setup times unnecessary, since additional table equipment is no longer required. All functions can be controlled with the foot switch. This allows a smooth transition between patient examination, education, and treatment. Numerous additional functions such as the ApexLocator or Implant Terminal in Sidexis 4 make Teneo the expert solution for every treatment situation.

Teneo supports all current reciprocating file systems, such as Wave One, Wave One Gold, Reciproc, and Reciproc blue.




Endodontics

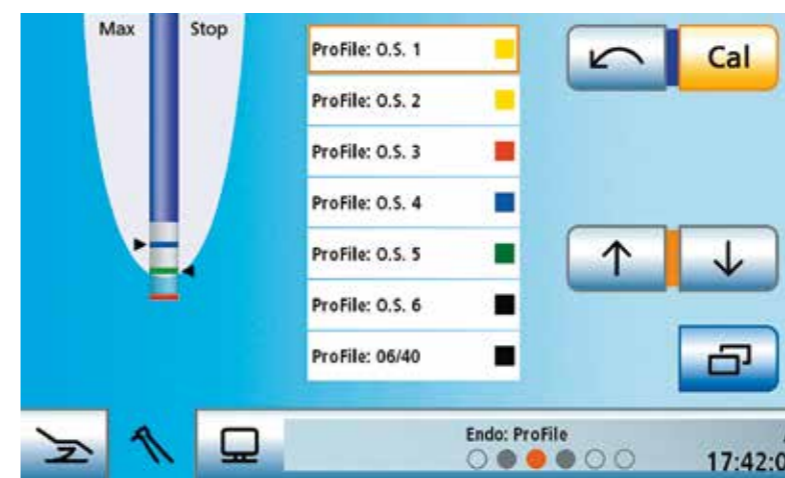

ApexLocator*

Precise torque control, the ApexLocator, and the unique file library streamline your workflow while offering maximum clinical safety. The only control you need is the Teneo foot switch.

Implantology

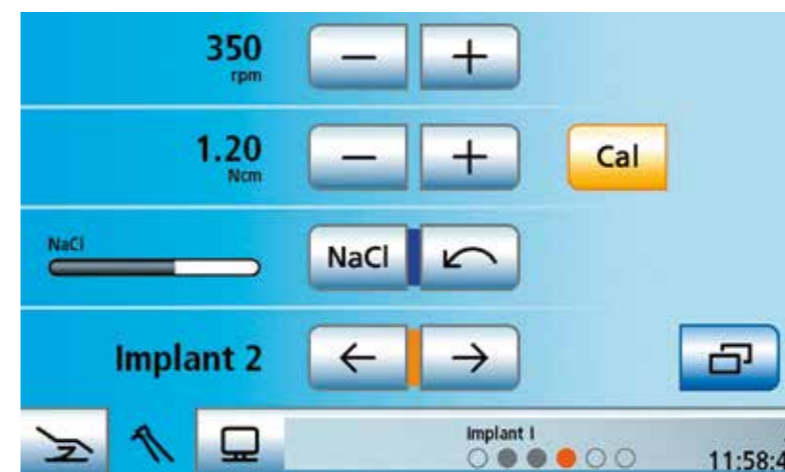

Implant Terminal*

Adjust Teneo to your specific needs manually or select an implant system from the integrated library. Important parameters such as speed, torque, and coolant flow are always under your control – efficient, sterile, and hands-free.

* ApexLocator and Implant Terminal are displayed on the Sivation monitor.

You want to be better in every detail – Teneo offers more

Teneo with its thoughtful details has everything you need to work perfectly in every situation. From a generously sized user interface to the flexible configuration of instrument positions, Teneo makes all the difference.



7"-EasyTouch display



The large touchscreen ensures intuitive and self-explanatory operation. It is customizable and upgradeable. In easy mode, it shows only those functions you need at any given time. Navigate safely and easily throughout the entire treatment.

6 + 1 instrument positions



Well-arranged and well-organized central tools and controls: Six instrument positions allow you to configure the treatment center, giving you even more ways adapt it to your requirements. The camera can be optionally placed on an auxiliary storage tray position

LEDview Plus operating light



LEDview Plus provides a perfect mixture of light intensity, light color, and light field size. With the integrated composite mode, light-curing filling materials can be processed at your own speed without premature hardening. The operating light can be optimally positioned and controlled hands-free via a sensor. Also available as ceiling light or as LEDview Plus/Heliodent Plus ceiling combination.

SiroCam AF+ camera



Optimize the visualization and documentation of your treatment: The SiroCam AF+ combines excellent image quality with easy operation. The image is displayed on the 22" monitor in high definition quality.



You create trust – Teneo helps you

Close communication with the patient is essential for first-class treatment results. A transparent treatment process contributes significantly to patient loyalty. Teneo allows clearly visualized on-site consultations, stores individual patient positions and thereby helps build trust in your operations.

Good advice from Teneo

- Use the integrated monitor to display intraoral images, 3D x-ray images or presentations using the patient communication menu.
- Create a trust-based relationship by explaining the treatment process to make the patient to understand the procedure and to reduce patient anxiety.
- Explain and document the treatment success with the SiroCam AF+ intraoral camera — of course with excellent high-quality images



Exceptional comfort

For many patients it is particularly important to feel appreciated and their individual needs are considered when it comes to dental treatments and the atmosphere in which it is performed. The comfortable technology of Teneo unobtrusively helps you achieve this. Patients who feel they are in good hands will gain confidence and stay faithful to you.

Why settle for less?

Whether old or young, big or small - with its need-based patient positioning, individual sitting positions that can be saved, massage function, lumbar support, slow-motion travel, and Thermo or Lounge upholstery, Teneo adapts to patient needs.



You know your patients – so does Teneo

A relaxed patient makes treatment easier for you and contributes to a perfect treatment result. Whether you require assisted seating for elderly patients or want to calm your patients' fear of dental treatment — Teneo exceeds the expectations of the most discriminating patient. Its massage and lumbar-support functions help your patients relax.

Needs-based posture

Positioning using the OrthoMotion function is particularly well suited for treating elderly or disabled patients: The backrest and seat retain the same angle, protecting the lumbar vertebrae.



Ergonomic patient positioning

The ErgoMotion function provides combined movements of the backrest and seat for an ergonomic extension of the spine. The patient is comfortably moved to the treatment position.



Shape and shade

To give your practice a unique atmosphere, Teneo provides a perfect combination of shape and shade, making it an absolute eye-catcher. Whether your preferred style is classic or modern, Teneo will always fit in elegantly with your office environment. This is also underscored by design awards such as the iF Design Gold Award or the nomination for the German Design Award.



Upholstery colors

The elegant world



Classic, timeless combinations for a truly sophisticated look.

The natural world



Soft, warm shades for a soothing treatment atmosphere.

The vital world



Vivid colors and bright tones for a modern, lively atmosphere.

* Available colors for Lounge upholstery.

Stitching colors



Standard Red Dark grey Light blue Rose

Paintwork colors



Uncompromising quality



Premium Thermo upholstery

Reduces the accumulated heat in the seat and back area and provides for patient relaxation as a result of this cooling effect.

Lounge upholstery

Its innovative material gives this upholstery option additional softness — for that bit of extra comfort.



Exquisite and functional

Give your practice a feeling of exclusivity. The Lounge upholstery provides a convincing synthesis of comfort and design. The soft material has a pleasant feel and is easy to clean and disinfect. Accentuated contrast stitching and the continuous footrest speak of high quality down to the smallest detail.

Teneo – safety made simple

The Teneo treatment center lets you meet the highest infection prevention standards – with smooth, easy-to-clean surfaces, easily removable functional components, and built-in sanitization adapters. Everything is geared to maximum infection control in minimum time with minimum effort to relieve your team of unnecessary work and to meet the highest hygienic standards.



Built-in sanitization adapter:

All instrument hoses can be connected to the integrated sanitization adapter for periodic cleaning of the waterways. The Auto-purge (daily rinsing) and Purge (rinsing after each treatment) programs flush the waterways automatically or disinfect them in the monthly sanitization cycles with the Dentosept water disinfectant. These processes are very safe and save time.



Easy-to-clean surfaces

Teneo has been designed from the ground up for efficient infection control management. All surfaces are smooth and can be cleaned and disinfected easily.



Thermodisinfestible parts

All infection control-critical parts are detachable and can be thermally disinfected, enhancing protection from cross-contamination.



Cuspidor

In the interest of ideal in-office infection control, the motorized cuspidor with its hands-free controls can be readily detached for easy cleaning.



Sanitization adapters for suction-hose cleaning

The suction hoses can be rinsed automatically at the press of a button, between treatments or even during lengthy treatment procedures. For maximum infection control, optional detergents can be added automatically.

Not just any Teneo – your Teneo

Teneo is more than a dental treatment center. It is an expression of your individuality. Whatever your favorite configuration, Teneo will always guarantee you one thing: interoperability in perfect harmony. Let yourself be inspired!

Teneo sets standards

The sophisticated ergonomics and functionality combined with quality materials help you mold your ideal personal workflow actively



Chair



Dentist element



Assistant element



Water unit



Sivision



Working stools

Your individual options

This is a complete listing of the standard equipment for the basic unit and all available individual options.

Chair	
Patient-specific chair positions*	■
Motorized headrest with automatic patient-specific headrest adjustment*	■
Premium Thermo upholstery	■
Lounge upholstery, by request with colored stitches	□
Lumbar support/massage function	■
Removable armrest, right armrest can be swiveled to side	□
Removable left armrest	□
Four-way foot switch	■
Wireless foot control	■
External device connection	□
Booster cushion for small patients	□
Adapter plate	□
OrthoMotion	■
ErgoMotion	■
SlowMotion travel	■

Dentist element (with EasyTouch, 6 positions, removable handles and instrument holder)	
Sprayvit M with light	□
1st drill drive (for motor or turbine)	■
2nd drill drive (for motor without turbine)	■
3rd drill drive (for motor without turbine)	■
4th drill drive (for motor without turbine)	□
BL electric motor	□
BL ISO C electric motor	□
BL Implant electric motor	□
Turbine hose	■
SiroSonic TL ultrasonic handpiece	□
Sirotom electrosurgery kit	□
Implantology function	□
Endodontic treatment function incl. reciprocating file function	□
ApexLocator with display on Sivision monitor	□
Integrated pump for saline solution	□
USB port	■
Hands-free operation with cursor control	■
Contra-angle handpiece T1 Line Endo 6 L with light	□
Contra-angle handpiece Implant 20:1	□

Assistant element (with 4 positions, removable holder)	
Satelec mini L.E.D. curing light	□
Surgical suction device	□
Sprayvit M with light	□
Suction device	■
Saliva ejector	■
Suction hoses with slide	□

Water unit	
Disinfection device with permanent disinfection and sanitization function	■
Suction hose cleaning with chemicals	■
Hydrocolloid port	□
Cuspidor valve for wet suction system	■
Automatic separator	□
Amalgam separator	□
Glass cuspidor, motorized swivel function, removable	■
Purge function with autopurge	■
Aqua function	■

Sivision digital	
SiroCam AF+	□
22" monitor on the tray, includes tray	□
22" monitor on the light support column	□
Camera connector in the dentist element	□
Cable set for PC connectivity	□

Working stools (long or short gas spring)	
Hugo Freehand (with foot release)	□
Hugo Manual (with manual release)	□
Hugo Manual Plus (with manual release and foot ring)	□
Carl Manual (with manual release)	□
Carl Manual Plus (with manual release and foot ring)	□
Saddle Stool Theo	□

Additional options	
Sirona tray	□
PC interface (Ethernet)	■
PC control with EasyTouch	■
Remote diagnostics	■
LEDview Plus operating light (triple-jointed support arm, with sensor)	□
Deckenleuchte LEDview Plus	□
LEDview Plus ceiling light	□
Adaptation for Heliodont Plus unit model	□

* Available with Sidexis 4.
 ■ In basic unit or as a package.
 □ Optional

Equipment options

Value packages



Innovation

- Sprayvit M for the dentist and assistant elements
- 2 BL electric motors (extra charge for BL ISO C)
- SiroSonic TL ultrasonic handpiece
- LEDview Plus operating light
- Up to two Hugo, Theo or Carl working stools at a discounted price



Endo

- Endodontic treatment function incl. reciprocating file function
- ApexLocator assembly
- T1 Line Endo 6 L contra-angle
- Basic Apex adapter (optional)
- Holder for ApexLocator



Implant

- Extended implantology function
- Integrated pump for saline solution
- BL Implant electric motor
- Implant contra-angle handpiece (20 : 1)



Sivision

- Monitor, 22"
- SiroCam AF+
- Monitor position on the tray or the light support column
- Camera connectivity in the dentist element



Booster cushion for small patients



Apex Locator



CEREC AI



External device connection



Integrated pump for saline solution



Implant adapter



Heliodont Plus model



22" monitor on the tray



22" monitor on the light support column

Dentsply Sirona for you

We support dentists in their quest for better, safer and faster dental care. If you routinely deliver maximum performance on a daily basis, you will expect no less from your equipment. This is why we provide innovative solutions that contribute to improving dental health. This means that you will always work to the highest standards of excellence with Dentsply Sirona. You can rely on that – today and in the future.

Hands-on innovation



As a global innovation leader for dental equipment, we continuously invest in research and in the future of modern dentistry. By networking digital technologies with integrated solutions and optimizing the treatment workflow, we create improved treatment results, more comfort and safety for the patient, and time and cost savings for the dental office.

What does this mean for Teneo?

Quality “Made in Germany”

The treatment center is developed, manufactured, and tested entirely at the Dentsply Sirona factory in Bensheim, Germany. The careful assembly on the premises justifies the highest expectations of Teneo, ensuring the best in long-term product quality and reliability.

Pioneering and innovative

Teneo is constantly evolving: The treatment center can be upgraded and has a USB interface for the integration of new functions and an Ethernet interface for easy troubleshooting.

High-quality workmanship

Only the best of materials are used for Teneo. The framework, chair and backrest support are made of a well-conceived combination of steel and aluminum. As a result, Teneo is lightweight yet extremely solid.

Procedural Solutions

Preventive
Restorative
Orthodontics
Endodontics
Implants
Prosthetics

Enabling Technologies

CAD/CAM
Imaging
Treatment Centers
Instruments

Dentsply Sirona

Sirona Dental Systems GmbH
Fabrikstraße 31, 64625 Bensheim, Germany
dentsplysirona.com

THE DENTAL
SOLUTIONS
COMPANY™

 **Dentsply
Sirona**