

GAC's Elation<sup>®</sup> MB  
Clear Polymeric-Metal Bracket

Strength, Beauty, Safety




# Our Commitment to Innovation: Elation® MB


Incorporated into the Elation® MB design is a durable, high-quality mechanical base. This innovative mechanical base has been designed to be reliably effective in both bonding and debonding procedures. In addition, this design improvement eliminates the need for a plastic primer.

The insert flows around the wire, which actually lowers friction, as the wire touches the insert only at four small points at the polymeric wall. Next, Elation® MB's steel translates the torque through the base without the usual distortion or give of plastic and other polymeric materials. Unlike inserts which only line the slot wall, Elation® MB's unique slot is placed within the body to strengthen the whole bracket, not just the slot. This placement also makes the steel insert harder to see, creating a more esthetic appearance.

We recommend our NeoBond® adhesive with Elation® MB; however, please feel free to utilize other adhesives. Elation® MB debonds cleanly and safely as well, ensuring your peace of mind. Elation® MB can be utilized with any arch wire, but for a total esthetic solution appliance, we recommend High Esthetic Wires.

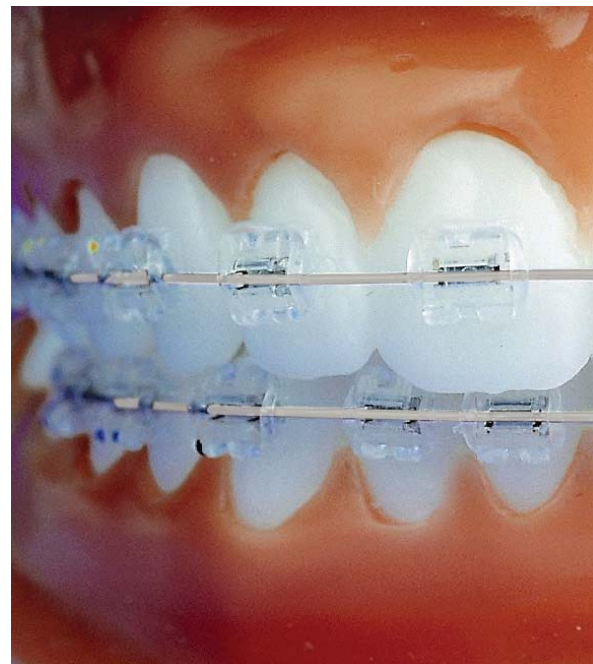

High Esthetic Wires shown

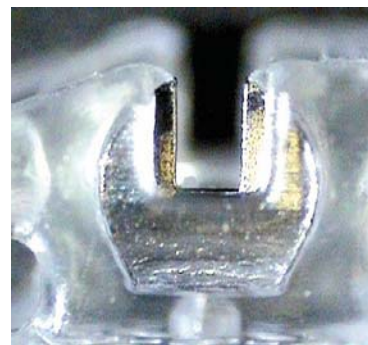

## BRACKETS

Maxillary Tooth	T	A	.018		.022	
			Right	Left	Right	Left
U1	12	5	133-111-60	133-211-60	133-112-60	133-212-60
U2	8	9	133-121-60	133-221-60	133-122-60	133-222-60
U3/hk	0	10	133-131-61	133-231-61	133-132-61	133-232-61
U4/5	-7	0	133-541-60	133-541-60	133-542-60	133-542-60
U4/5 hk	-7	0	133-541-61	133-541-61	133-542-61	133-542-61

Mandibular Tooth	T	A	Right	Left	Right	Left
L1-2	0	0	133-611-60	133-611-60	133-612-60	133-612-60
L3/hk	-11	7	133-331-61	133-431-61	133-332-61	133-432-61
L4	-17	0	133-641-60	133-641-60	133-642-60	133-642-60
L4/hk	-17	0	133-641-61	133-641-61	133-642-61	133-642-61
L5	-22	0	133-651-60	133-651-60	133-652-60	133-652-60
L5/hk	-22	0	133-651-61	133-651-61	133-652-61	133-652-61

KITS (brackets only)	.018"	.022"
UL5-5 CS & BC HK	K133-551-62	K133-552-62

CONSULTATION MODEL	Reference Number
Elation MB	TYPI33-777-77

NeoBond Adhesive Kit	52-800-00
NeoBond Paste Syringe, 4g	52-800-10
Gel Etch Syringe w/ tips, 5ml (37% phosphoric acid)	52-800-16
Primer, 6ml	52-800-11

To learn more, visit us online at  
[dentsplysirona.com](http://dentsplysirona.com).



Dentsply GAC International  
One CA Plaza, Suite 100  
Islandia, NY 11749

©2017 Dentsply GAC Intl.  
All Rights Reserved  
0280-MM-120-133-01 Rev. 02