




Low Profile Buccal Tubes

Providing a smooth, comfortable patient experience.

Low Profile Buccal Tubes

LP buccal tubes with integrated hooks are metal injection molded from high quality medical grade stainless steel providing a smooth, comfortable patient experience. This popular design comes in a full assortment for first and second molars.

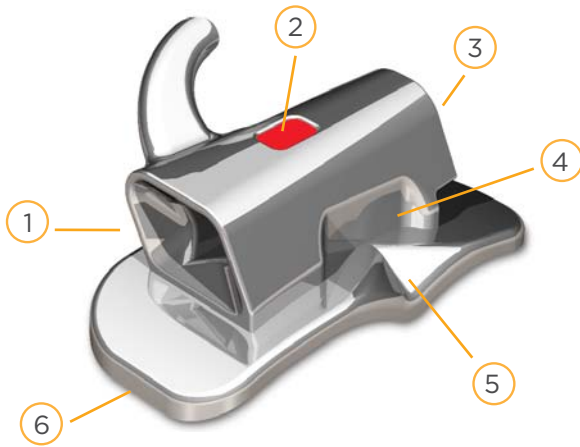


- 1. Enhanced opening allows for easier wire insertion
- 2. Color recess for simple ID
- 3. Expanded exit port for reduced friction and wire impingement
- 4. Instrument notch for a secure grip during placement
- 5. Occlusal indent on base pad offers easier placement in the buccal groove.
- 6. 80 gauge foil mesh base

Low Profile mini Buccal Tubes

LP mini Tubes are specifically designed for partially erupted 2nd molars. The mini design provides a solution for early stages of bonding. The bonding base features a gingival offset with a mesially-extended contour ideal for adapting to a limited bonding surface.

- 1. Enhanced opening allows for easier wire insertion
- 2. Low profile hook designed to enhance patient comfort and minimize cheek irritation
- 3. Color recess for simple ID
- 4. Two-piece mesh construction designed to provide consistent and predictable bonding
- 5. Gingivally offset with a mesially-extended base to facilitate placement on partially erupted second molars

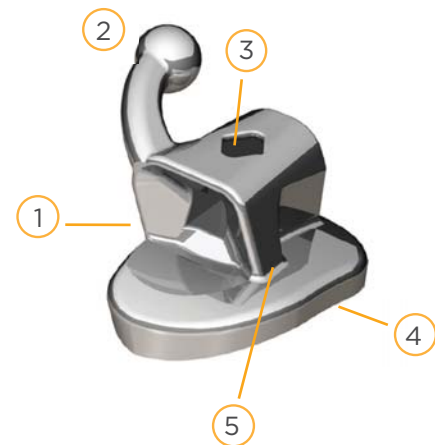

Dentsply GAC International
One CA Plaza, Suite 100
Islandia, NY 11749

©2017 Dentsply Sirona
All Rights Reserved
0280-MM-120-168-02 Rev. 05