



OmniArch® PLUS by GAC

Precision Driven
Performance





Remastering a Masterpiece

When it comes to performance, precision and predictability go hand in hand. Manufactured at our state-of-the-art facility in Sarasota, Florida, USA, OmniArch PLUS takes everything you love about the OmniArch system—and takes it to the next level. Reimagined and rebuilt from the base up, OmniArch PLUS pushes the limits of precision to deliver the ultimate in performance.

Precision PLUS Performance

- Small enough for the esthetic conscious patient
- Large enough for efficient usability
- Integrated ball hook for stability and malleability
- Triple chamfered slot for easy wire placement
- 80 gauge mesh base for predictable bonding and debonding

OmniArch[®] PLUS

OmniArch[®] PLUS Advantage



Torque In Base

Compound Contoured Base

Smooth Swept Tie-Wings



MIM (Metal Injection Molded) Design


Delivers 50% increase in accuracy of torque and 33% in offset for greater consistency, accuracy and precision.



Center Scribe Line

Provides a visual reference for alignment with the long axis of the clinical crown.

Rhomboid Shape

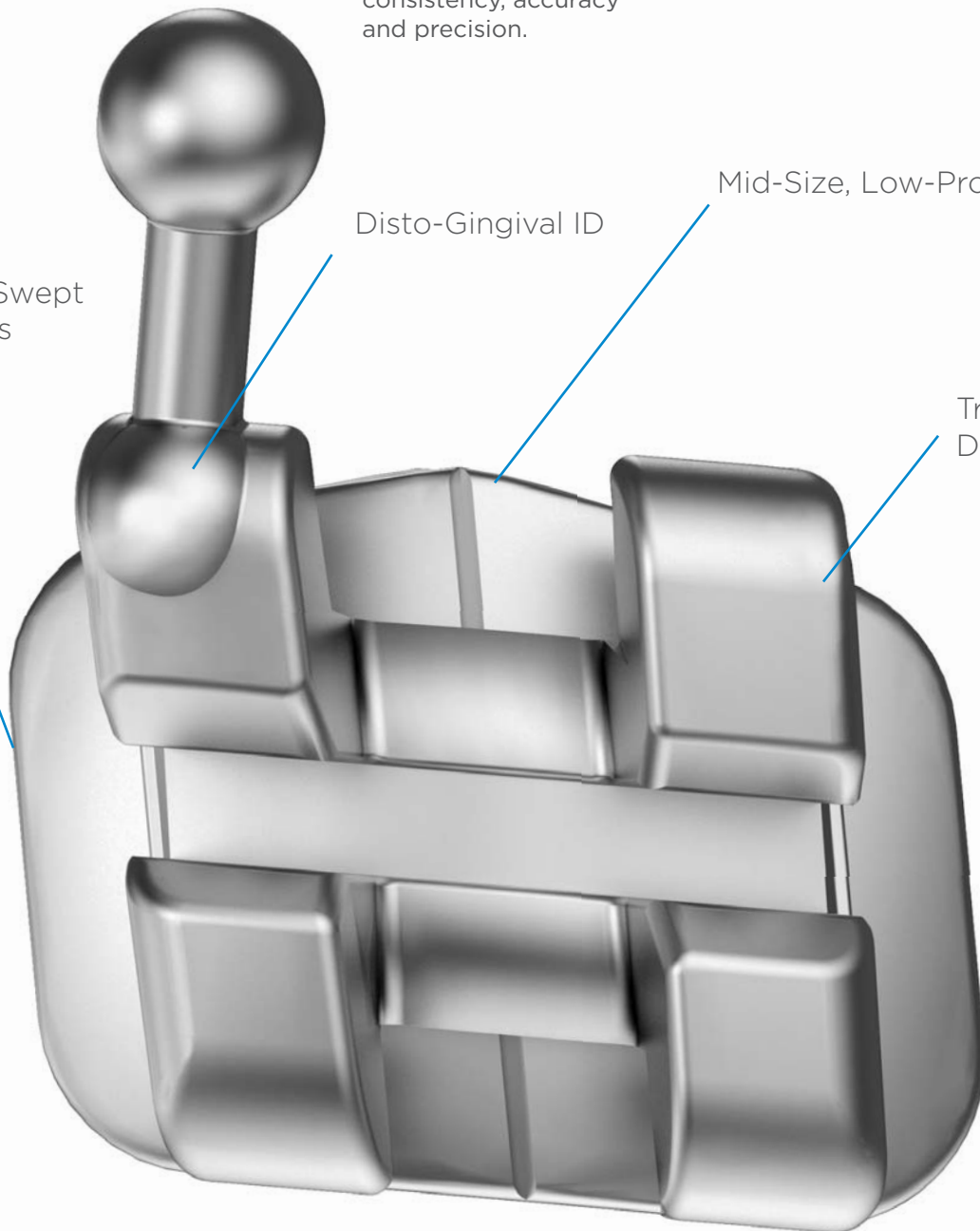

Disto-Gingival ID

Mid-Size, Low-Profile Twin

True Twin Design



80 Gauge Mesh Design

80 gauge single mesh, compound contoured base engineered for a precise anatomical fit.



Integrated Ball Hook Design

Enhanced design consistency while maintaining the essential functionality, such as malleability and strength.



Triple Chamfered Slot Walls

Easy engagement, designed to reduce the chance of the archwire binding or crimping.

Engineering Success


It's all-but-impossible to improve upon an industry standard like the Omni without a good idea of what you want to accomplish. Every design enhancement to OmniArch PLUS was subjected to an exhaustive explorative process. Each idea was tested, refined and then tested again. Nothing was considered complete until it met or exceeded our exacting standards.

Insight PLUS Innovation

- Metal injection molded for the ultimate in production consistency
- Fine-tuned offset and torque for enhanced control
- Unique polishing process for smooth, comfortable finish
- Welded base body for stability and precision
- Digitally engineered tie-wings for outstanding workability

(Available in Base Rx*, Euro (*MBT) Rx and Vertical Slot Rx.)

*Dentsply Sirona Orthodontics Base Rx has prescription values that are equivalent to the Roth prescription. ROTH is a registered trademark of Roth Licensing, LLC. All rights reserved. MBT is a registered trademark of 3M Unitek.

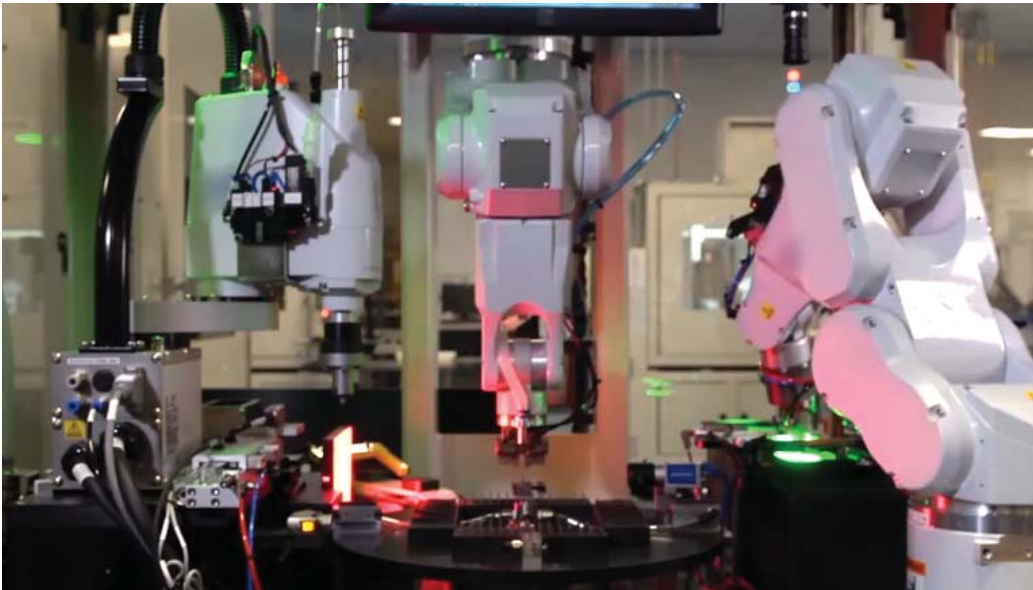

The Future of Orthodontics is Here

The OmniArch PLUS is the first complete prescription to be produced in Dentsply Sirona's state-of-the-art manufacturing facility in Sarasota FL. The brackets are metal injection molded, a manufacturing process that allows for incredibly intricate design elements to be reproduced with a high degree of accuracy. This means that all of the design elements and bracket programming are engineered to be as close to perfect as possible.

OmniArch® PLUS by GAC



Schedule Your Free Lunch and Learn with GAC!

We'll bring lunch to your office and provide a hands-on demonstration for you and your staff!¹

To sign up for your FREE GAC Lunch and Learn, visit:

dentsplysirona.com/orthodontics/ortholearn

1-800-645-5530

¹ Some geographical restrictions may apply.

The Dental Solutions Company

Dentsply Sirona is the world's largest manufacturer of professional dental products and technologies. We create industry-leading positions and platforms across consumables, equipment, technology, and specialty products. We are committed to introducing new, innovative products and complementary, end-to-end solutions to drive better, safer and faster dentistry.

To learn more, visit us online at dentsplysirona.com.



Dentsply GAC International
One CA Plaza, Suite 100
Islandia, NY 11749

©2017 Dentsply Sirona
All Rights Reserved
0280-MM-120-J77-01 Rev. 03